



## 2026 MD LASER IPL + RF THERAPY TREATMENT

Skin rejuvenation

Pigment therapy

Vascular therapy

Acne treatment





# THREE TREATMENT MODES

The three treatment modes to meet a variety of treatment needs.




# INTELLIGENT OPERATING SYSTEM

The operator only needs to select the part and skin color to be treated, and the system will automatically generate the recommended treatment parameters. Even an inexperienced operator can easily grasp the operation of the instrument



# SHR "In-motion" Hair Removal

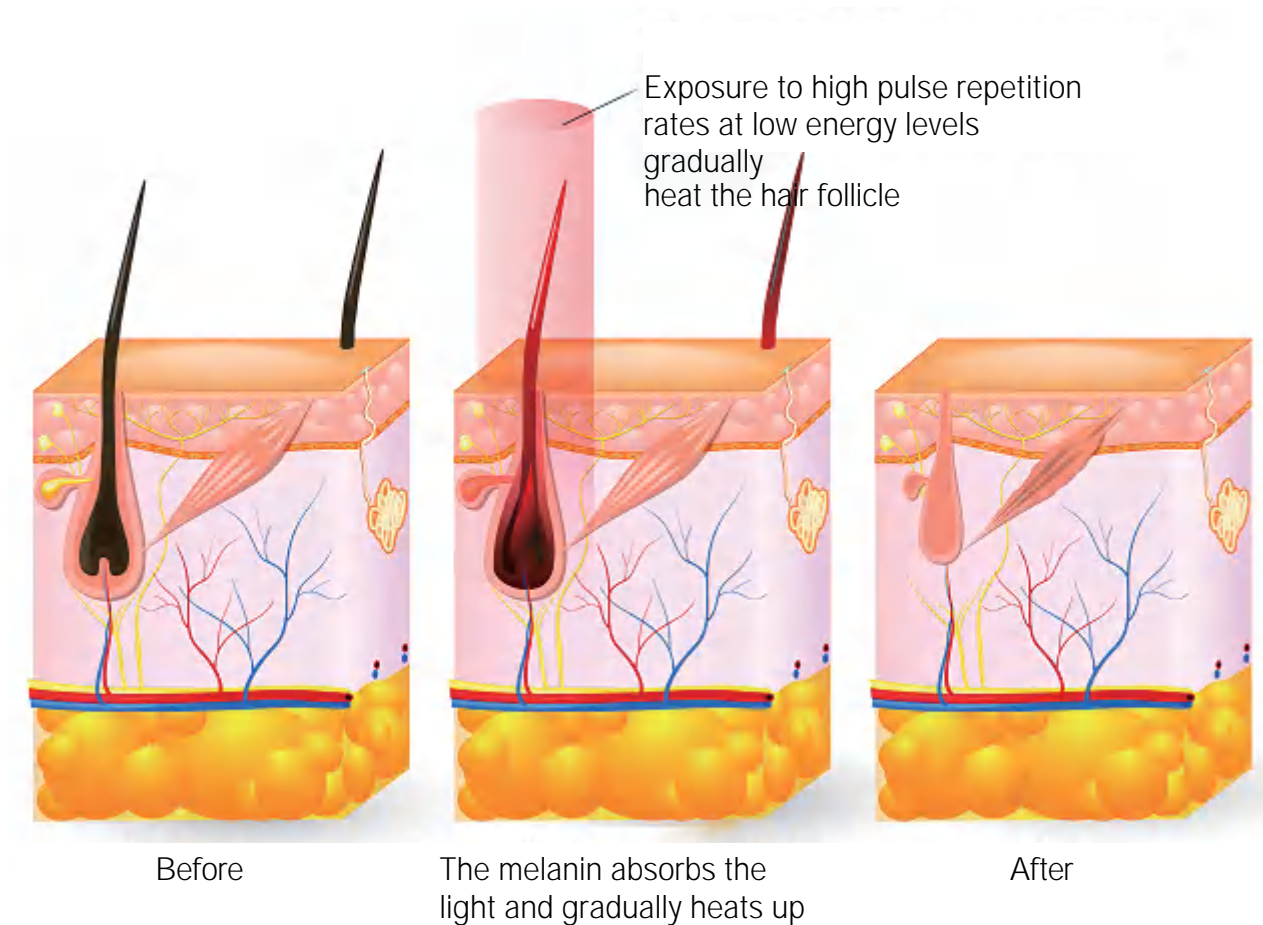

SHR stands for Super Hair Removal and it's the latest breakthrough innovation in IPL technology to launch onto the ever-advancing beauty and aesthetics market.



## Faster

SHR in-motion treatment method has much improved the treatment speed.

The SHR machine repetition rate is up to 10Hz. The handpiece could shoot 10 shots in 1 second.

## More Comfortable


Unlike traditional IPL treatments, which use high energy levels that may cause discomfort and pain, SHR Skin Pulse treatments are mild and comfortable, SHR scans the treatment zone, gradually heating the skin with repetitive low fluence pulses, ensuring that you'll experience comfort even in the most sensitive of areas.

## All Skin Types

The gradual heat distribution of SHR means treatments are totally comfortable, non-irritating and virtually pain free making it ideal for even sensitive dark skin types.




IPL SKIN TREATMENT

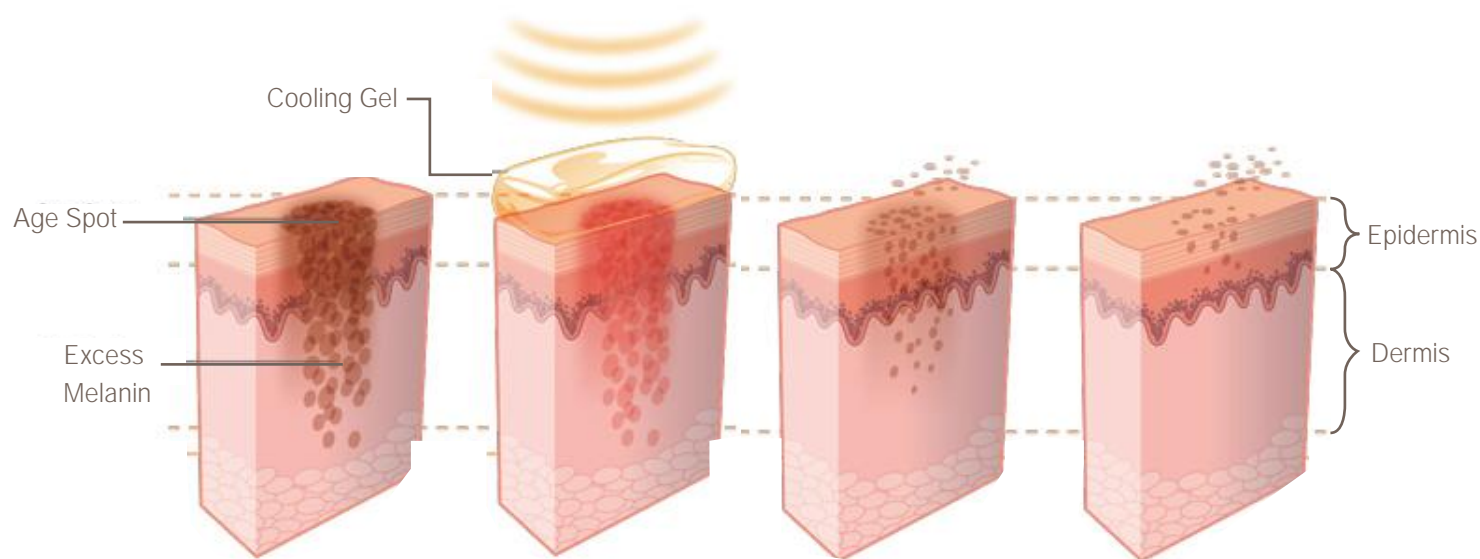

IPL technology is a type of Intense Pulsed Light that generates intense bursts of red or blue light.

These heat the dermal tissue without harming or removing the outer epidermal layers of the skin. The body's natural process then work to break down and remove these marks.

The appearance of any skin imperfection is then visibly reduced.

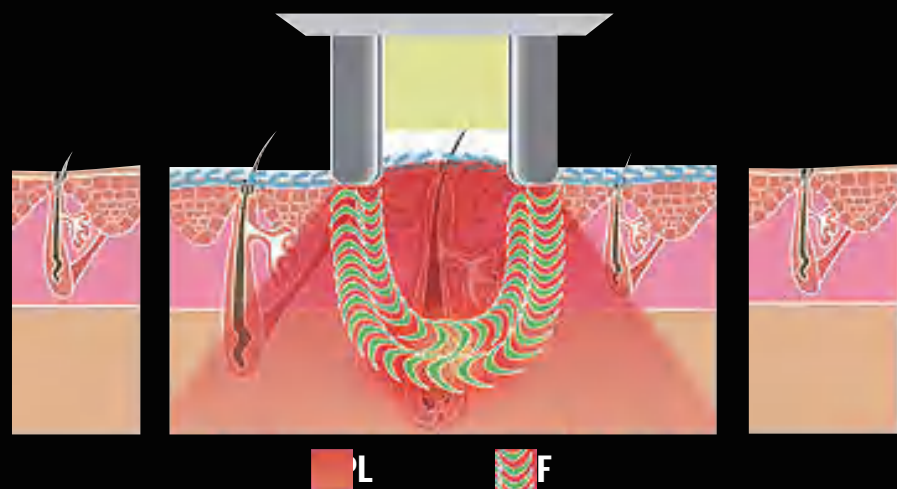


The IPL+RF, while doing IPL treatment, radio frequency energy is output from the metal splints on both sides of the IPL treatment head, which has the effect of lifting and tightening.




## REAL-TIME INTELLIGENT MONITORING SYSTEM

20 intelligent detection functions, including water flow, water temperature, handle temperature, handle connection, water level etc. Ensuring the stability of the machine and find out the possible problems of the machine in time

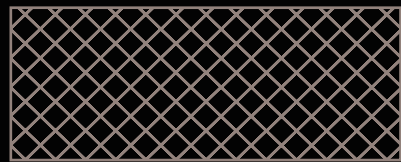

**15\*50mm larger spot size**

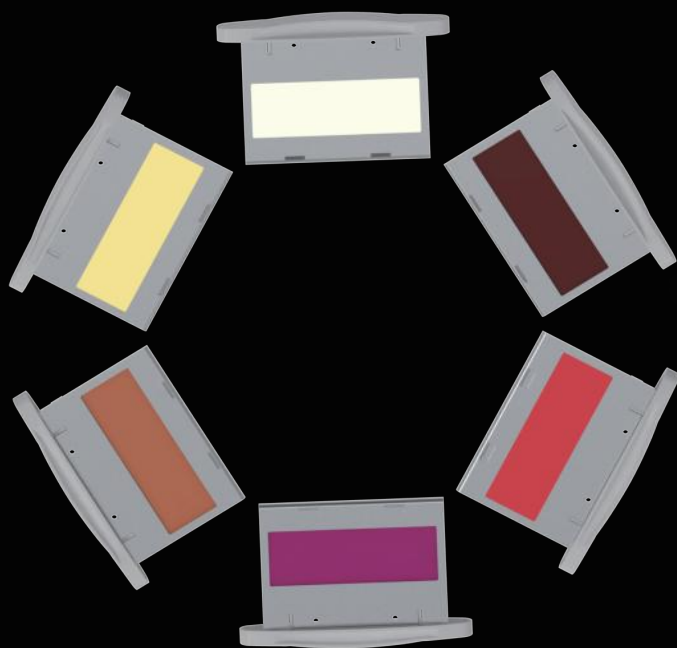

Cooling sheet work detection, etc

handle connection status detection

Water level detection

Water level detection

### 2000W High power

The output energy is high enough to ensure that the energy will not decay even in the fast light mode, ensuring a good treatment effect

handle connection status detection

### Strong Cooling System

Ensuring handle surface is always cool and the treatment always comfortable.


At the same time, making sure that the machine works continuously for 24

hours without damage by the high temperature



# Six different filters

One handle is equipped with 6 filters with different wavelengths, which can meet various treatment needs.

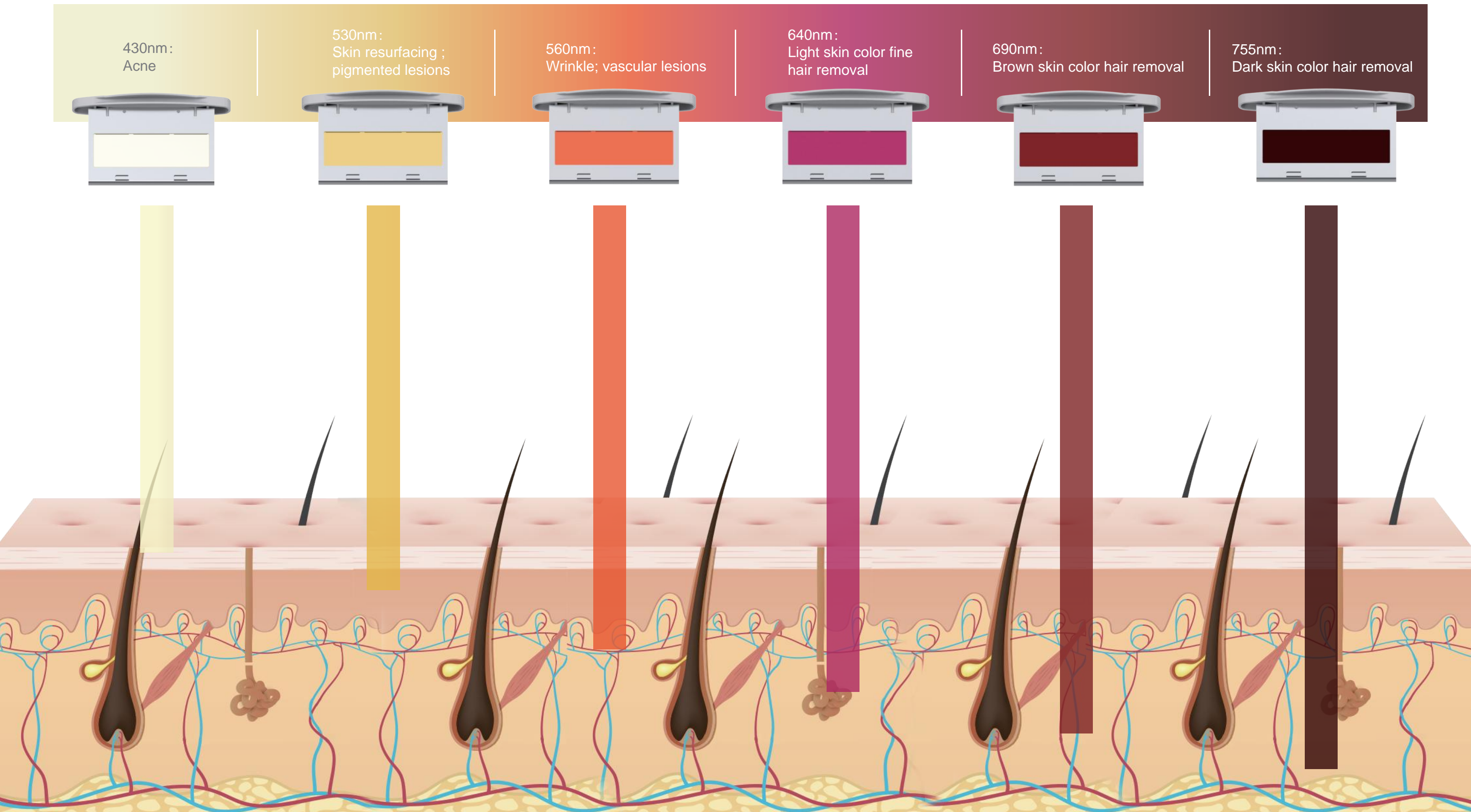

## Handle Cooling

The handle system adopts 108W high-power refrigeration film, which is twice as much as other machines on the market, ensuring comfort-able and safe treatment


## Silver Reflector

The reflector is made of silver, which is not easy to be oxidized, has higher energy and longer life


## Sapphire Crystal



The sapphires crystal integrated into the handpiece that contacts the skin during the laser/IPL treatment. It cools and protects the epidermis and reduces pain and discomfort. It also allows better light transmission than metal contact cooling.

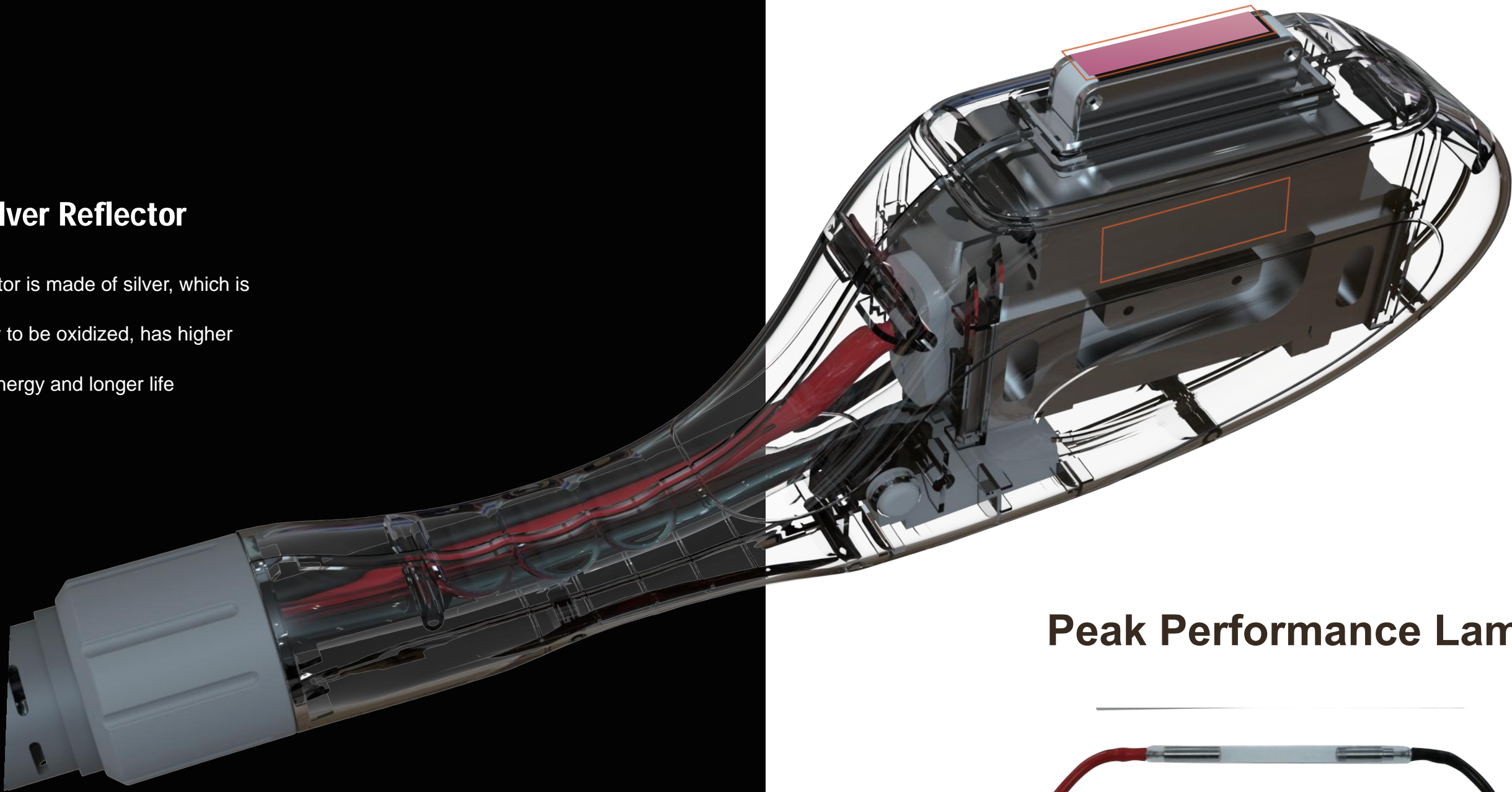

## Peak Performance Lamp

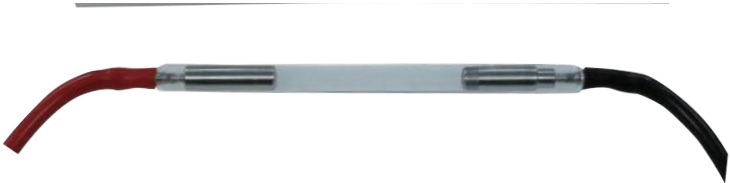



# SPECIFICATIONS

Power	1500W
Wavelength	430nm,530nm,560nm,640nm,690nm,755nm
working modes	IPL, SHR, E-light
IPL energy	10-60j/cm2
Spot size: SR	15*50mm
SHR energy	3-15j/cm2
SHR frequency	1-10Hz
SHR pulse width	1-10ms
RF power	160W
RF energy	1-50j
Dimensions	23wX16dX16h 100LBS
Voltage	AC220V±10% 10A 50HZ , 110v±10% 10A 60HZ

